

# AMBU® TRAINING MANIKINS

Features	AmbuMan School	Ambu Uniman+	AmbuMan Basic	AmbuMan Airway Instrument	AmbuMan Airway Wireless	AmbuMan Advanced (+ IV)
Mouth to Mouth ventilation	•	•	•			
Mouth to nose ventilation	•	•	•			
Bag-valve-mask ventilation		•	•	•	•	•
Use of Guedel airways				•	•	•
Airway management: Et tube, laryngeal tube & mask and other SGAs				•	•	•
Jaw-thrust manoeuvre			•	•	•	•
Gastric insufflation				•	•	•
Ambu hygiene system	•	•	•			
Insufflation bag (Ambu hygiene system)	20	25	25			
Manikin face mask (Ambu hygiene system)	•	•	•			
Chest compression	•	•	•	•	•	•
Max compression depth (cm)	6	7.5	7.5	7	7	7
Max insufflation - litres	0.4	1.8	1.2	1.2	1.2	1.2
Hand position display		•	•	•	•	•
Infinitely adjustable chest stiffness				•	•	•
Pulse				manual	manual	electrical
Placement of defibrillation electrodes	•	•	•	•	•	•
Defibrillation						•
I.V. access - need IV trainer 34041				optional	optional	optional
Wireless CPR software					•	•
ECG display						•
Carrying bag	carton box	•	•	•	•	•
Weight (kg)	0.8	4.5	6.0	13	13	13.3
Multilingual manual-languages*	17	17	12	10	10	19

\*All manikins have multilingual manual in: GB, FR, IT, ES, PT, DE, DK, HU, PL. 34040, 34045, 34046, have also: CZ, GR, NO, SE, SI, TR, JP, CN. 34042, 34048 have also NL. 34045 has also FI.

## • 34045 AMBU® MAN SCHOOL

Training manikin for education in school of all procedures for basic life support resuscitation. It can be used for realistic training of external chest compressions and the placement of AED pads. The unique hygienic system prevents cross contamination when demonstrating and training of the principals of ventilation through mouth and nose.

### Guideline conform


It allows to perform CPR training following the latest international CPR guidelines.

### Compact and portable

It has a low weight and compact packaging.

### Easy AED training

Due to the closed shaped torso, defibrillation pads can be placed to perform a realistic AED training.

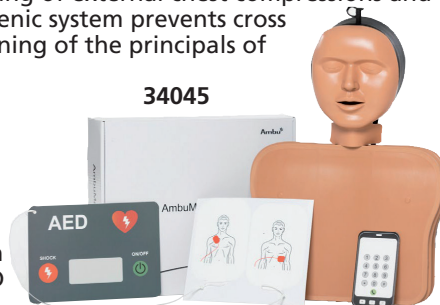

## • 34046 AMBU® UNIMAN+

Inexpensive training manikin with all the well-known and sought-after Ambu look-listen-feel features. The closed shape of the torso with its anatomical correct landmarks enables students to learn the complete algorithm of Basic Life support, including the use of Automatic External Defibrillators (AED) and the placement of the AED Pads.

### Open Airway


The airway is opened only when the head is tilted - there are realistic airway resistance and chest rise.

### Hygienic System

Unique patented hygienic solution as the Ambu Man with disposable head bags and washable face pieces.

### Compact Bag

The head of the Ambu Uniman+ can be placed inside the manikin. Everything is packed in a small bag.

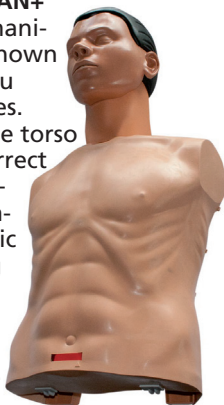


34046

## • 34047 AMBU® MAN BASIC

Cost efficient and reliable training manikin designed to fulfill all requirements for the BLS algorithm training. The closed shaped torso with its anatomical correct landmarks allows trainees to learn and train the use of Automatic External Defibrillators (AED) and the placement of pads. It allows an optimal training of chest compression and ventilation with realistic resistance and expiration through mouth and nose. The unique hygienic system prevents cross contamination. The correct compression depth and required ventilation volume can be monitored by the feedback indicator on the manikin.

### Ventilation as in real life

The trainees will feel the realistic airway resistance and the chest will start rising when the head is correctly tilted and air blown into the manikin. Size: 70x37 cm.

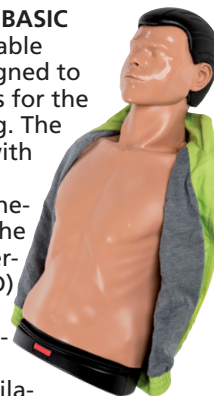

34047

## • 34039 HEAD BAG - box of 100 pcs

## • 34040 AMBU® MAN ADVANCED

Professional training manikin made to fulfil all the needs of Advanced Life Support (ALS) training. The focus of the manikin is on the algorithm: Airway Management, ECG simulation, defibrillation and the training of i.v. access technique. Integrated Ambu Wireless Technology: the browser-based Ambu Manikin Management Module allows the wireless control, documentation and monitoring of the trainee's performance through the central interface, where all data and information come together: compression depth, correct hand position, ventilation volume, stomach inflation.

### Customised scenario training

The intuitive user interface allows a quick and easy creation of custom scenarios and adapt them to different course concepts, in accordance to the demand of instructors and the current ERC / AHA guidelines. The CPR parameters and manual events are documented wirelessly and included into the evaluation and debriefing.

### Individual trainee evaluation in group

More efficient training, each participant receives his personal results (up to four participants).

- 34041 AMBU IV TRAINER
- 34029 BATTERY PACK

## • 34042 AMBU® MAN AIRWAY I

## • 34048 AMBU® MAN AIRWAY W with CPR

Adult-sized training manikins that provide an optimal lifelike representation of the human anatomy. Instant feedback on ventilation volume, stomach inflation, chest compression depth in mm and wrong hand position.

### PC connection (34048 only)

of Ambu CPR Software Kit, Ambu CPR Printer or USB Interface.

### Make intubation possible

with ET tube, Guedel, Laryngeal Mask and combi-tube.


Ambu Man Airway includes torso with hooks for arms, computer interface with CPR software kit (34048 only), track suit jacket, training mat and carrying bag. CPR software kit selectable languages: GB, FR, IT, ES, DE, PT, PL, DK, HU, NL, RU, JP.

- 34043 COMPLETE SET: ARMS, LEGS, TROUSERS AND BAG
- 34041 AMBU IV TRAINER
- 34029 BATTERY PACK



ERC  
Guideline  
conform  
2015

34040  
34041

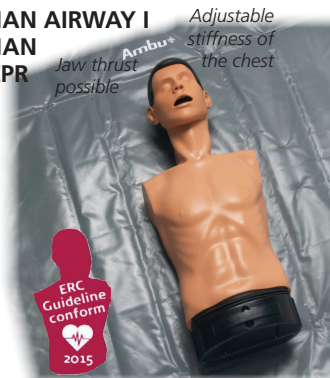

ERC  
Guideline  
conform  
2015

34042,  
34048

Size: 80x35 cm  
Weight: 13 kg



34042,  
34048